



دستگاه قالب گیری بافت (بلوکه کردن بافت) تیشو امبدینگ مدل TE100

- دارای مخزن ذوب پارافین از جنس استیل ضد زنگ با ظرفیت ۳/۵ لیتر
- توانایی تنظیم دمای مخزن ذوب از دمای محیط تا +۱۰۰ درجه سانتیگراد
- توانایی تنظیم صفحه گرم کننده از دمای محیط تا +۱۰۰ درجه سانتیگراد
- دارای هات پلیت از جنس آلومینیوم آنادایز شده ابعاد ۳۲ × ۲۷ سانتیمتر
- دارای نازل خروجی مجهز به شیر برقی قابل کنترل با فوت سوئیچ
- دارای فوت سوئیچ (پدال پای) برای استفاده آزادانه از دو دست
- دارای کشتیو زیر صفحه گرم کننده برای جمع آوری پارافین اضافی
- دارای المنت گرمکن (هیتر) جداگانه برای شیر نازل خروجی پارافین
- دارای هلدر فورسپس شش خانه از جنس آلومینیوم و چراغ روشنایی LED
- دارای کنترل پنل برای تنظیم دیجیتال صفحه گرم کننده و مخزن ذوب پارافین
- دارای نمایشگر LED برای مشاهده دمای تنظیم شده و دمای واقعی
- دارای حافظه داخلی برای ذخیره کلیه تنظیمات دستگاه
- دارای یکسال گارانتی و ده سال خدمات پس از فروش

TE100 Tissue Embedding Unit


TE100 is composed by paraffin container, dispensing unit, heated working area. The ergonomically elevated heated working area is illuminated with LED light and allows comfortable processing of specimen.

TE100 LCD display of control panel; you can visualize all the temperatures. The user is allowed to set the parameters through the menu function.

TE100 easy-to-operate control panel consists of a membrane keyboard with push button for instrument

TE100 is equipped with:

- Paraffin melting tank made of stainless steel, max capacity 3.5 liter
 - Hot plate to perform embedding operations, equipped with canals to collect and drain the paraffin in excess to drawer, made of aluminum
 - Paraffin dispenser with metering knob for a constant flow rate regulation.
- The working area is properly illuminated by a cold light source
- Foot switch for paraffin flow activation without using operator's hands
 - Heating forceps holders on both sides of Paraffin dispenser

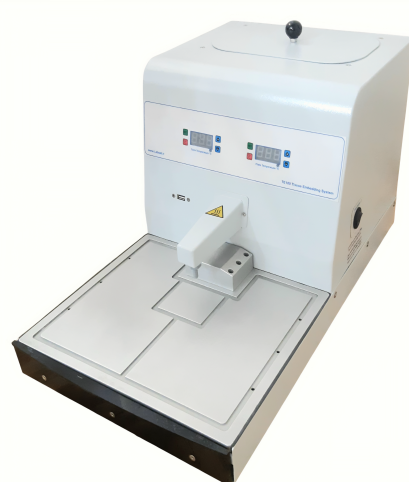

تلفن: ۰۲۱- ۸۸۳۲۷۷۳۰ فکس: ۰۲۱- ۸۸۴۹۴۲۶۶ تلفن همراه: ۰۹۱۲-۵۳۷۳۲۷۸

www.labset.ir sales@labset.ir

@Labset.ir کانال تلگرام و صفحه اینستاگرام